

## DIGITAL THERMOMETER FOR CHILDREN

SC

44 flex



### Digital Thermometer SC 44 flex - Penguin


- flexible measuring tip
- temperature reading in 10 secs
- last temperature value record
- fever alarm
- beep
- waterproof.

The SC 44 flex is a fast-maximum thermometer, i.e. the temperature measurement is displayed after approx. 10 seconds. The final measurement is indicated by a beep signal.

The SC 44 flex has an integrated fever alarm. At temperatures above 37.5 ° C, the thermometer is producing short signal tones in succession and the LCD display flashes.

With memory function, the last measured value is saved and auto displayed after power on.

The thermometer comes within a plastic case inside the product packaging box.

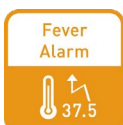

EAN 4022356016036  
Art.-Nr.: 01603

More info at [www.scala-electronic.de](http://www.scala-electronic.de)

CE  
0197

## DIGITAL THERMOMETER FOR CHILDREN

SC

44 flex

### Specifications

Type	Maximum Thermometer
Measuring range	32.0°C to 43.9°C Temperature below 32.0°C: display "L" for low. Temperature above 43.9°C: display "H" for high.
Measuring duration	appr. 10 seconds
Resolution	0.1°C
Deviation	± 0.1°C between 34°C and 42.0°C at 18°C to 28°C ambient temperature. ± 0.2°C between 32.0°C and 33°C and between 42.0°C and 43.9°C at 18°C to 28°C ambient temperature.
Fever alarm	Above 37.5°C (display is flashing)
Display	LCD-Display, 3 digits
Shutdown	appr. 3 minutes after the end of measurement
Storage temperature	Temperatures: -10°C to 55°C Relative humidity: 15% to 95%
Battery info	1 x LR41 (1.5V), capacity 30mAh, weight 0.6g
Battery life	min 36 months
Battery indicator	yes
Water resistant	yes
Colour	white with motif
Dimension	154 x 21 x 12 mm
Device weight	appr. 18 gr
Country of origin	China
Operational manual languages	German, English, French, Dutch, Czech, Polish, Italian, Danish, Swedish, Spanish



### Logistics details

	Quantity in pcs / kg	Nett weight (Kg)	Gross weight (Kg)	Dimensions (cm)
Single package	1	0.022	0.052	19x2.5x6 cm
Sales unit	6	0.132	0.365	16.5x6.8x23 cm
	120	2.6	8.1	36x50x36 cm
One layer pallet	960	22	73.2	120x80x54 cm
Whole pallet	3840	88	292.8	120x80x200 cm

EAN 4022356016036

Art.-Nr.: 01603

More info at [www.scala-electronic.de](http://www.scala-electronic.de)

CE  
0197